MANUAL DE COLHEITAS DE AMOSTRAS

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Collection precautions:

Donor should NOT eat, drink, smoke or chew gum for 30 minutes before collecting saliva sample.

Ensure the sponge tip does NOT come into contact with any surface prior to collection.

Intended use: ORAc collect•Dx is intended for use in the non-invasive collection of saliva samples. Human DNA from the saliva sample is isolated, stabilized, and suitable for use in FDA cleared molecular diagnostic applications. Saliva samples collected using ORAc collect•Dx are stabilized and can be transported and/or stored long-term at ambient conditions.

Contents: Kit contains stabilizing liquid.


Warnings and precautions:

- Rx only. For professional use.
- Read all instructions carefully prior to sample collection; deviation may result in inadequate sample and impact DNA yield.
- Choking hazard. Caution should be used when inserting sponge into the mouth.
- Donors with xerostomia (dry mouth) may not collect adequate sample using these instructions resulting in lower DNA yield and an invalid sample for use.
- Wash with water if stabilizing liquid comes in contact with eyes or skin. Do NOT ingest. See MSDS at www.dnagenotek.com


Storage: 15°C-25°C

Label legend:


	Catalog number
	In vitro diagnostic medical device
	CE marking
	UKCA marking
	Manufacturer
	Medical device
	Do not re-use
	Consult instructions for use
	Collect sample by (Use by)
	Lot number
	Authorized Representative

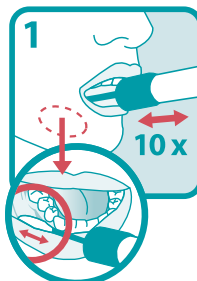
Open package and remove collector without touching sponge tip. Place sponge as far back in the mouth as comfortable and rub along the lower gums (see close-up image) in a back and forth motion. Gently rub the gums 10 times. If possible, avoid rubbing the teeth.

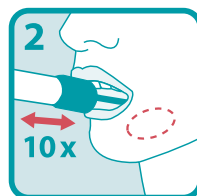
Gently repeat rubbing motion on the opposite side of the mouth along the lower gums for an additional 10 times.

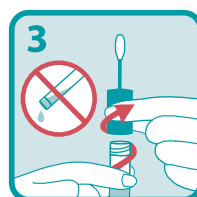Hold the tube upright to prevent the stabilizing liquid inside the tube from spilling. Unscrew the blue cap from the collection tube without touching the sponge.

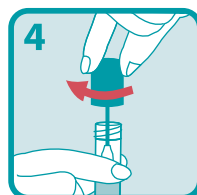
Turn the cap upside down, insert the sponge into the tube and close cap tightly.

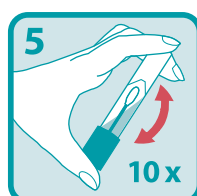Invert the capped tube and shake vigorously 10 times.

USER INSTRUCTIONS

Read all instructions prior to collection

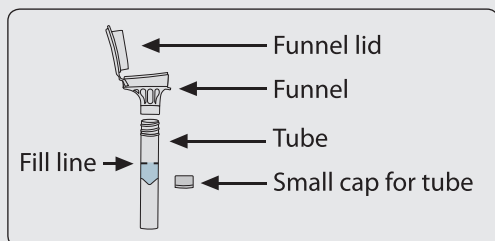Collection precautions:

Do NOT eat, drink, smoke or chew gum for 30 minutes before giving your saliva sample.

Do NOT remove the plastic film from the funnel lid.

Intended use: For the collection of human DNA from saliva samples.

Contents: Kit contains stabilizing liquid.

Warnings and precautions: Wash with water if stabilizing liquid comes in contact with eyes or skin. Do NOT ingest. See MSDS at www.dnagenotek.com.

Small cap may pose a choking hazard.

Storage: 15°C / 30°C

Summary and explanation of the kit:


Oragene•DNA is a self-collection kit that provides the materials and instructions for collecting and stabilizing saliva specimens.

Label legend:


	Consult package insert
	Collect saliva by (Use by)
	In vitro diagnostic medical device
	Catalog number
	CE Marking
	Caution, consult instructions for use
	Storage instructions
	Authorized Representative
	Manufacturer
	Lot number

Procedure:


Most people take between 2 and 5 minutes to deliver a saliva sample following steps 1 to 5.

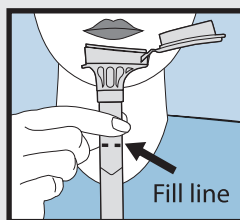
- 1 Spit into funnel until the amount of liquid saliva (not bubbles) reaches the fill line shown in picture #1.

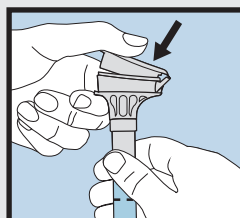
- 2 Hold the tube upright with one hand. Close the funnel lid with the other hand (as shown) by firmly pushing the lid until you hear a loud click. The liquid in the lid will be released into the tube to mix with the saliva. Make sure that the lid is closed tightly.

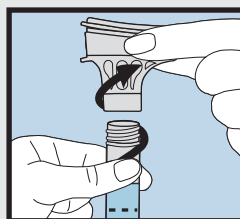- 3 Hold the tube upright. Unscrew the funnel from the tube.

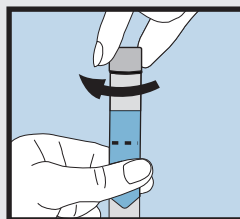- 4 Use the small cap to close the tube tightly.

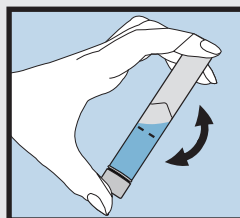
- 5 Shake the capped tube for 5 seconds. Discard or recycle the funnel.

NOTE: The sample must be collected immediately after opening the transparent bag.

NOTE: Do not force or excessively press or bend the swab when taking the sample from the patient, as this may result in accidental breakage of the shaft. The swab shaft has a moulded breakpoint in proximity of the tip, which allows intentionally breaking the swab inside a test tube.

1. Put on gloves, open the bag and remove the plastic tube (Figure 1A).
2. To facilitate the opening of the tube cap, hold the swab by the plastic tube with one hand and with your thumb press the top of the white cap (Figure 1B).
3. Holding the swab by the grip (white cap), remove the swab from the plastic tube and place the tube on a clean surface (Figure 1C).
4. Collect the sample. During sampling, the swab end (tip) shall only come in contact with the collection site, so as to reduce contamination risks.
5. Once you have taken the sample, put the swab back into the plastic tube and securely close it (Figure 1D).
6. Identify the collected sample. Send the sample to the laboratory for testing as soon as possible.

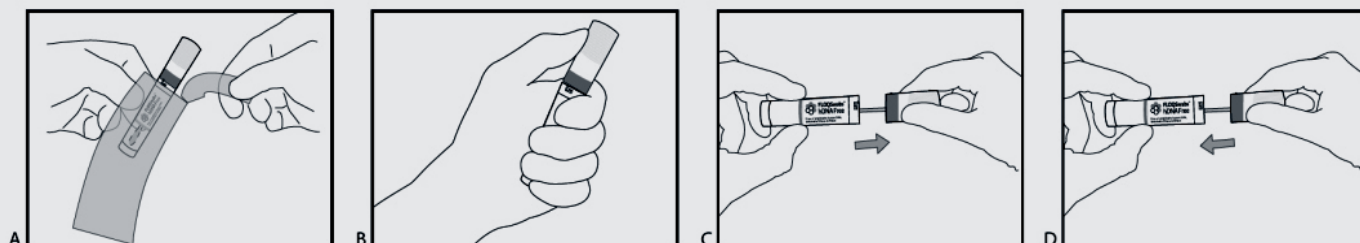WASTE DISPOSAL

Unused swabs may be considered as non hazardous waste and disposed of accordingly.

Dispose of used swabs as well as any other contaminated disposable materials following procedures for infectious or potentially infectious products. It is the responsibility of each laboratory to handle waste and effluents produced according to their nature and degree of hazardousness and to treat and dispose of them (or have them treated and disposed of) in accordance with applicable regulations.

Table of Symbols

Symbol	Meaning
	Manufacturer
	Do not reuse
	Catalogue number
	Temperature limits
	Use before

Symbol	Meaning
	Consult the instructions for use
	Batch code (lot)
	Contents sufficient for <n> tests
	Do not use if the packaging is damaged
	Non-sterile

Copyright © 2019 Copan Italia S.p.A.
All rights reserved

 Copan Italia S.p.A.
Via F. Perotti, 10
25125 Brescia - Italy
Tel +39 030 2687211
Fax +39 030 2687250

Email: info@copangroup.com
Website: www.copangroup.com

North American Distributor:
Copan Diagnostics Inc.
26055 Jefferson Avenue
Murrieta, CA 92562 USA
Tel: 951-696-6957
Fax: 951-600-1832

 **COPAN**
Innovating Together™